

Automated biochemical analyzer

SPOTCHEM™ EZ

SP-4430 | Operating Manual

Chapter 6

APPENDIX


6-1 Transmission Specifications

6-1-1 Protocol

6-1-2 Format

6-1-3 Format for measurement result (format 1)

6-1-4 Format for measurement result (format 2)

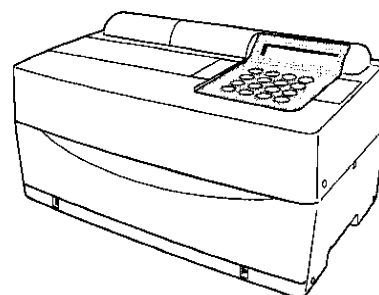

6-1 Transmission Specifications

Chapter 6 APPENDIX

6-1-1 Protocol

6

Transmission format	Start-stop system (asynchronous) unidirectional transmission by serial transmission format (in compliance with JIS X5101).
Data format	<p>One character consists of the following 11 bits.</p> <p>Start bit: 1 bit</p> <p>Data bit: 7 bits (ASCII code)</p> <p>Parity bit: 1 bit (even)</p> <p>Stop bit: 2 bits</p>
Baud rate	<p>Baud rate can be selected out of the following 6 rates.</p> <p>300, 600, 1200, 2400, 4800, 9600bps</p>
Hand shake	Suppression by CFT and RTS are possible. (Suppression is not set by default.) XON/XOFF control is not available.
Time gap	2-second waiting time is inserted between each block (from < ETX > or < ETB > to < STX >).
Forced finish	Data transmission is sometimes forced to stop halfway by key operation. It is not promptly stopped by pressing a key, but transmission continues until < ETX > or < ETB > is output.

6-1-2 Format

Block structure is regular. One block consists of start, data and end. This is illustrated below in the following descriptions.

Start	Data	End
--------------	-------------	------------

- **Start (S)**

Start of each block is < STX >.

Start of block is expressed as S in the following illustrations.

- **Data**

Data (text) of each block is the main body of transmission contents, and is expressed by arrangement of ASCII characters.

< CR >, < LF >, < RS > or < US > is sometimes involved in data. Characters other than these cannot be involved.

- **End (E)**

End of each block is < RTX > or < ETB >.

< ETX > or < ETB > is distinguished by whether it is in the last block or not.

If it is in the last block, it is < ETX >, and if it is in the middle block, it is < ETB >.

The block end is expressed by E in the following illustrations.

6-1-3 Format for measurement results (format 1)

The measurement result (format 1) is the same as the "regular format" in SP-4410 or SP-4420. The receiving program designed to receive the measurement results of SP-4420 (regular format) can normally receive the measurement results (format 1) of SP-4430.

■ Transmission data of measurement results (format 1)

When transmitting the measurement results with format 1, the number of blocks differs depending on the combination of the reagent strips.

A. When only Multi Reagent Strips are measured.

S	Header	Multi Reagent Strip measurement results	E
---	--------	---	---

B. When only Single Reagent Strips are measured.

S	Header	Single Reagent Strip measurement results	E
---	--------	--	---

C. When Multi and Single Reagent Strips are measured.

S	Header	Multi Reagent Strip measurement results	E	S	Single Reagent Strip measurement results	E
---	--------	---	---	---	--	---

● Format of header

001	002	003	004	005	006	007	008	009	010	011	012	013	014	015	016	017	018	019	020	021	022
■	■	/	■	■	/	■	■						■	■	:	■	■			CR	LF
023	024	025	026	027	028	029	030	031	032	033	034	035	036	037	038	039	040	041	042	043	044
I	D	#		■	■	■	■	■	■	■	■	■	■							CR	LF

Start	Finish	Description
001	001	Measurement date. Year (The last two digits of year), month (1~12), date (11~31). No zero control. YMD format is always applied regardless of the date setting.
014	018	Measurement date. Time (01~23), Minute (01~59). No zero control.
027	036	When ID is available, ID is output. When there is no ID, only the first 10 digits are output. The measurement number is expressed in 4 digits without zero control. For the measurement number, blank is used for 0311~036.

●Format of Multi Reagent Strip measurement results

001	002	003	004	005	006	007	008	009	010	011	012	013	014	015	016	017	018	019	020	021	022
	M	U	L	T	I	:	■	■	■	■	■	■	■	■	■	■				CR	LF
023	024	025	026	027	028	029	030	031	032	033	034	035	036	037	038	039	040	041	042	043	044
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
045	046	047	048	049	050	051	052	053	054	055	056	057	058	059	060	061	062	063	064	065	066
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
067	068	069	070	071	072	073	074	075	076	077	078	079	080	081	082	083	084	085	086	087	088
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
089	090	091	092	093	094	095	096	097	098	099	100	101	102	103	104	105	106	107	108	109	110
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
111	112	113	114	115	116	117	118	119	120	121	122	123	124	125	126	127	128	129	130	131	132
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
133	134	135	136	137	138	139	140	141	142	143	144	145	146	147	148	149	150	151	152	153	154
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF

Start	Finish	Description
008	017	Name of Multi Reagent Strip
023	027	Item name
029	029	Abnormal mark ● Within the range of standard values: Blank (20H) ● Low value: <US>(1FH) ● High value: <RS>(1EH)
030	034	Measurement value
036	041	Unit symbol
042	042	Temperature ● 37℃: Blank(20H) ● 30℃: " + " ● 25℃: " * " ● Items other than enzyme: Blank
045	154	The same repetition as 023-044. When the number of items is less than 6, the redundancy makes up for the blank (20H).

●Format of Single Reagent Strip measurement results

001	002	003	004	005	006	007	008	009	010	011	012	013	014	015	016	017	018	019	020	021	022
■	■	■	■	■	■	■														CR	LF
023	024	025	026	027	028	029	030	031	032	033	034	035	036	037	038	039	040	041	042	043	044
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
045	046	047	048	049	050	051	052	053	054	055	056	057	058	059	060	061	062	063	064	065	066
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF
067	068	069	070	071	072	073	074	075	076	077	078	079	080	081	082	083	084	085	086	087	088
■	■	■	■	■		■	■	■	■	■	■		■	■	■	■	■	■	■	CR	LF

Start	Finish	Description
001	007	When only Single Reagent Strips are measured: Blank (20H) When Multi Reagent Strips are measured: Fixed character string "SINGLE"
023	027	Item name
029	029	Abnormal mark ● Within the range of standard values: Blank (20H) ● Low value: <US>(1FH) ● High value: <RS>(1EH)
030	034	Measurement value
036	041	Unit symbol
042	042	Temperature ● 37℃: Blank(20H) ● 30℃: " + " ● 25℃: " * " ● Items other than enzyme: Blank
045	048	The same repetition as 023~044. No extra output is made. The block length of this block changes according to the number of Reagent Strips (items). For example, if the number of the Reagent Strips are 2 (2 items), the block length is completed as 66 bytes.

●Data errors

In the Multi Reagent Strip and Single Reagent Strip measurement results, data parts for each item (22 bytes including end CR/LF) are output as follows when measurement errors occur.

●Error in range or prozone (OVER)

023	024	025	026	027	028	029	030	031	032	033	034	035	036	037	038	039	040	041	042	043	044
■	■	■	■	■		■	■	■	■	■	■	■		■	■	■	■	■	■	CR	LF

Start	Finish	Description
023	027	Item name
029	035	<ul style="list-style-type: none"> ● Over the range: Fixed character string "OVER >" ● Under the range: Fixed character string "UNDER >" ● Prozone (OVER): Fixed character string "OVER >"
037	041	<ul style="list-style-type: none"> ● Over the range: Upper limit of measurement range ● Under the range: Lower limit of measurement range ● Prozone (OVER): Upper limit of measurement range
042	042	Temperature <ul style="list-style-type: none"> ● 37°C: Blank(20H) ● 30°C: "+" ● 25°C: "*" ● Items other than enzyme: Blank

●Prozone (CANT' MEAS) or calibration error by magnetic card

023	024	025	026	027	028	029	030	031	032	033	034	035	036	037	038	039	040	041	042	043	044
■	■	■	■	■		■	■	■	■	■	■	■	■	■	■	■	■	■	■	CR	LF

Start	Finish	Description
023	027	Item name
029	041	Error message <ul style="list-style-type: none"> ● Prozone (CANT' MEAS): "CANT' MEAS " ● Magnetic card, CAL error, L: "CAL. ERROR L1" ● Magnetic card, CAL error, H: "CAL. ERROR H1"
042	042	Temperature <ul style="list-style-type: none"> ● 37°C: Blank(20H) ● 30°C: "+" ● 25°C: "*" ● Items other than enzyme: Blank

6-1-4 Format of measurement results (format 2)


The measurement result (format 2) is the same as the "extended format" in SP-4410 or SP-4420.

The receiving program designed to receive the measurement results of SP-4420 (extended format) can normally receive the measurement results of SP-4430 (format 2).

■ Transmission data of measurement results (format 2)

A. When only Multi Reagent Strips are measured.

When transmitting the measurement results with format 2, 1 item is output as 1 block. The order of the items are, Multi Reagent Strip, Single Reagent Strip.

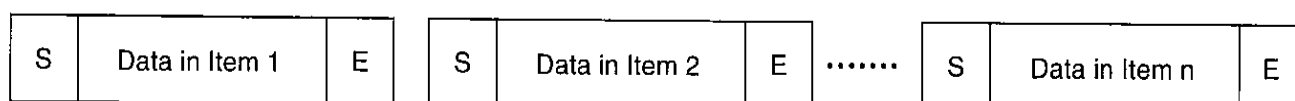

※ The end of the block is all block <ETX>.

●Format of "data in item x"

001	002	003	004	005	006	007	008	009	010	011	012	013	014	015	016	017	018	019
■	■	■	■	■	■	■	■	■	■		■		■	■	■	■	■	
020	021	022	023	024	025	026	027	028	029	030	031	032	033	034	035	036	037	038
■	■	■	■	■	■	■	■		■	■	■	■	■		■		■	CR
040	041	042	043	044	045	046	047	048	049	050								
■	■	■	■	■	■	■	■	■	■									
051	052	053	054	055	056	057	058	059	060	061	062	063	064	065	066	067	068	069
■	■	■	■	■		■	■	■	■	■		■	■	■	■	■	■	
070	071	072	072	074	075	076	077	078	079	080	081	082	083	084	085	086	087	088
0																		CR

Start	Finish	Description
001	010	When ID is available, ID is output. When there is no ID, only the first 10 digits are output. The measurement number is expressed in 4 digits without zero control. For the measurement number, blank is used for 031~036
012	012	Sample number types. For ID, "1". For number, "0"
014	015	The maximum number of items in the measurement
017	018	The order of the maximum number of items
020	027	Measurement date. Year (The last two digits of year), month (1~12), date (1~31) No zero control. YMD format is always applied regardless of the date setting
029	033	Measurement date, Time (0~23), Minute (0~59). No zero control
035	035	Data error. Normal=0, Under the range=1, Over the range or prozone error (OVER)=2, Low value=3, High value=4, Prozone error (CAN'T MEAS)=5, Calibration error by magnetic card L1=6, Calibration error of magnetic card H1=7
037	037	<ul style="list-style-type: none"> ● 37°C: "0" ● 30°C: "2" ● 25°C: "1" ● Items other than enzyme: "0"
040	049	Name of Multi Reagent Strip (For Single Reagent Strip, blank (20H) is used)
051	055	Item name
057	061	Measurement value. When the data error=1, lower value, the data error=2,5,6,7, upper values are applied
063	068	Unit symbol